




# One Step Alcohol Saliva Test




Wondfo One Step Alcohol Saliva Test Strip is a rapid, sensitive test to detect the presence of alcohol in saliva for blood alcohol concentration (BAC) greater than 0.02%. The test requires no instrumentation, calibration or special training to be used effectively.

The Alcohol Saliva Test Strip is ideal for Zero Tolerance Testing. It is suitable for testing employees at workplaces and agencies with zero tolerance policies!

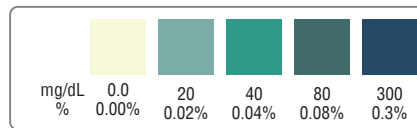

Color Chart

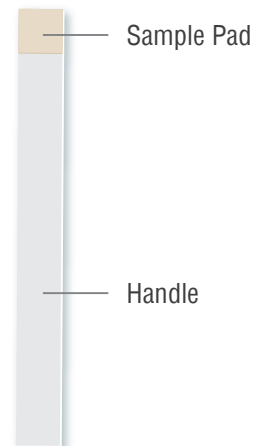


- Easy to use
- 18 month shelf life
- Results in only 2 minutes
- Approximate blood alcohol concentration

## Test Procedure



- 1** Take the strip test out of the sealed pouch.



- 2** Add one drop of saliva directly to the reaction pad.

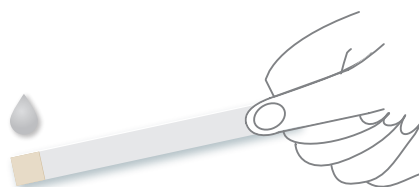

- 3** Read the result at 2 minutes by matching the color on the test strip to the color chart.

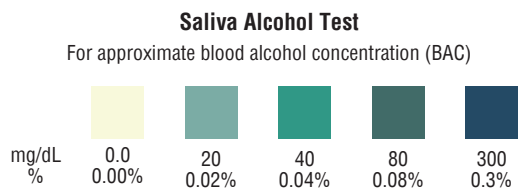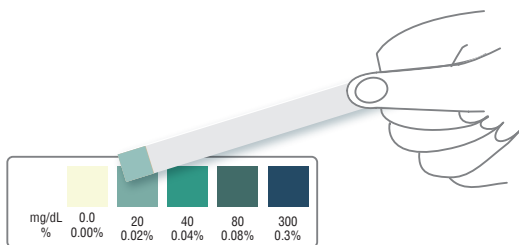

## Ordering Information

Item Code	Description	Kit Size / Case Size	Complexity
W53-S-100BX	100-Strips; Alcohol Saliva Test Box 100 Pcs	100 / 25	CLIA WAIVED
W53-S-25BX	25-Strips; Alcohol Saliva Test Box 25 Pcs	25 / 168	CLIA WAIVED